

SAM 2PX

The SAM 2 PX articulator has all of the basic features and functions of the SAM 2P articulator but with additional enhancements. A rigid centric locking mechanism allows the instrument unrestricted sagittal movement opening. In addition, condylar inclination and Bennett guidance angulation are adjustable. A wide variety of options and accessories are available, including different condylar curved pathways (remo-

vable inserts) and Bennett guidance curves (removable inserts) both of which allow the instrument to be fully adjustable. The protrusive adjustment is built-in. Transparent condylar housing covers keep the upper and lower members of the articulator together when the centric lock is disengaged and also allows for eccentric movements.



SAM 2PX Features



Spring loaded cover centric locking device with secondary security screw



Bennett angulation indicator and setting mechanism




Condylar Housing inclination angle indicator



Incisal pin with numeric scale



Facebow earpiece pin

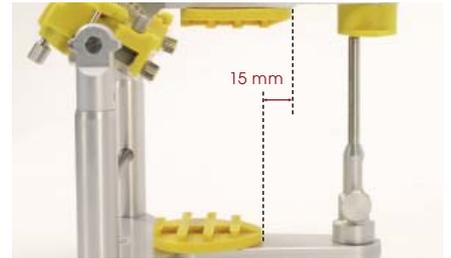

Mandibular mounting plate 15 mm sagittal posterior position



Selection of cast mounting plates, MPS magnetic or screw type. (AXIOSPLIT® optional)



Vertical Support Rod for upper member of articulator.




Tilt support rods for angular positioning of articulator



Protrusion adjustment up to 5.5 mm. Retrusion on lateral side 15 degrees



Options: Micro Adjustable Incisal Pin

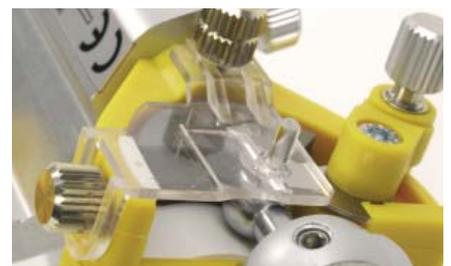

Transparent condylar housing cover



Options: Adjustable Incisal Table



Options: Bennett and Condylar Insert Set PX



Visible high quality and precision assures reliability and longevity